



One drop, to unify.

The safe and effective solution for wound healing and skin lesions haemostasis.



EFFECTIVE ¹⁻¹³

- Rapid wound closure (60-90 sec¹) ¹⁻⁸
- Excellent aesthetic results ^{2-6,9,10}
- Infection risk prevention ^{2,3,8,10-12}
- Trauma ^{2,10,11,13} and inflammation ^{7,11} reduction



EASY ^{1,2,5,10,11,15}

- Ready to use
- Easy to apply



PRECISE ¹⁴

- Product fluidity,
- Ergonomic bottle,
- Dedicated tip,

allow for a controlled application even in areas that are difficult to treat¹⁴ especially for closing wounds on the face and the neck or where the application is difficult (e.g. curved areas, such as elbow, knee, etc.).



VERSATILE ¹⁴

- 3 sizes: 0.25 ml, 0.35 ml, 0.50 ml



COMFORTABLE ^{2-4,10-13}

- Painless ^{2,7,11,15}
- Transparent and water resistant ^{5,10,15}
- No dressing ^{1,2,4,10,13}
- Spontaneous detachment 5-8 days after application ^{4,5,7,10,11,13}
- Indicated in children and in uncooperative patients ^{2,6,10}



COST-EFFECTIVE BENEFITS ^{2,4,5,10,16}

- No product waste
- Follow-up visits not necessary ^{1,3-5,7,10,12,13}
- No anaesthesia ^{2,5,10,15}



0.25 ml


0.35 ml

0.50 ml



**HAEMOSTASIS ALSO ENSURED
IN COAGULOPATHIC AND
ANTICOAGULATED PATIENTS** ⁸


The advantages of two monomeric adhesives in one drop.

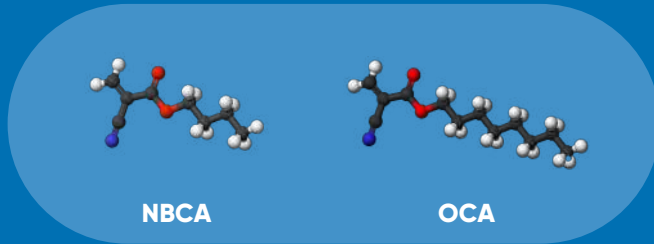


(N-BUTYL-2 CYANOACRYLATE) (2-OCTYL CYANOACRYLATE)

The two monomers (NBCA+OCA) ensure:

- > Elasticity ^{8,10}
- > No polymerisation in the presence of air ¹⁴
- > Polymerisation temperature of 45 °C, much lower than monomeric skin adhesives ¹⁴

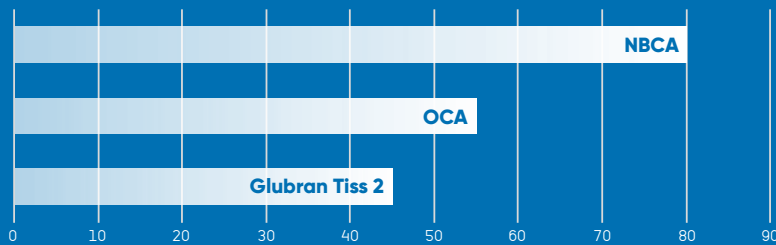

- > High tensile strength (27.34 N) ¹⁴

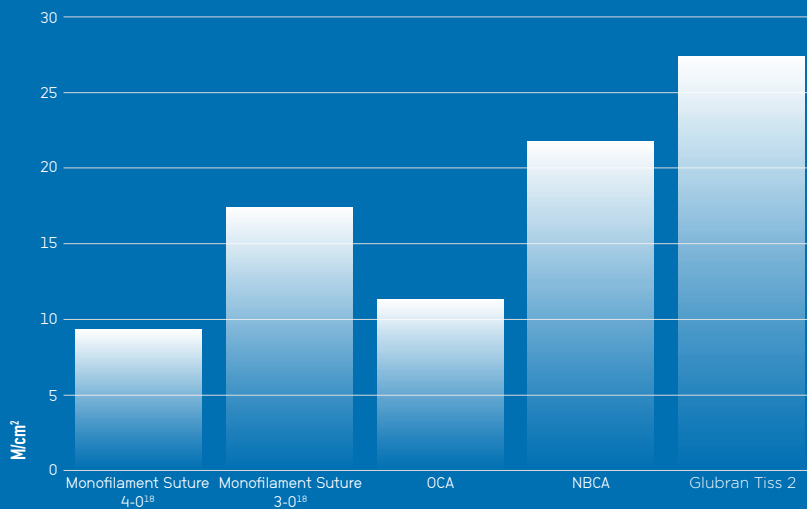

AUTHORISED FOR PAEDIATRIC USE ^{16,17}

GLUBRAN® Tiss 2

Treatment of lacerations, wounds and skin injuries

**on the face, limbs and
torso caused by trauma
or surgical incisions.**



Emergency Room

SKIN LESIONS CAUSED BY TRAUMA AND CUTS

Trauma laceration



Before



After 4 weeks

Operating Room

IATROGENIC WOUNDS

C-Section



Cutting in Thyroid Surgery



Trocar insertion



Certified Quality.

Glubran® Tiss 2 has obtained the CE Product Certification in 2007 from the 0373 Notified Body – Istituto Superiore di Sanità (Italian Health Institute).



GEM s.r.l. has obtained the following Quality Certifications from the TÜV SUD Certification Body:



EN ISO 13485:2016

UNI EN ISO 9001:2015

Medical Device Single Audit

Program (MDSAP)

Bibliography

1. Lombardi CP, Bracaglia R, Revelli L, Insalaco C, Pennestrì F, Bellantone R, Raffaelli M. Aesthetic result of thyroidectomy: evaluation of different kinds of skin suture. *Ann Ital Chir.* 2011 Nov-Dec;82(6):449-55
2. Gulalp B, Seyhan T, GURSOY S, Altinors MN. Emergency wounds treated with cyanoacrylate and long-term results in pediatrics: a series of cases; what are the advantages and boards? *BMC Res Notes.* 2009 Jul 14;2:132.
3. Chibbaro S, Tacconi L. Use of skin glue versus traditional wound closure methods in brain surgery: A prospective, randomized, controlled study. *J Clin Neurosci.* 2009 Apr;16(4):535-9.
4. Gennari R, Rotmensz N, Ballardini B, Scevola S, Perego E, Zanini V, Costa A. A prospective, randomized, controlled clinical trial of tissue adhesive (2-octylcyanoacrylate) versus standard wound closure in breast surgery. *Surgery.* 2004 Sep;136(3):593-9.
5. Graefe H, Wollenberg B, Brocks C. Application of skin adhesives in head and neck surgery: analysis of cosmetic results, applicability and costeffectiveness of cyanoacrylate-based adhesives. *HNO.* 2008 Sep;56(9):947-54.
6. Craven NM, Telfer NR. An open study of tissue adhesive in full-thickness skin grafting. *J Am Acad Dermatol.* 1999 Apr;40(4):607-11.
7. Farion K, Osmond MH, Hartling L, Russell K, Klassen T, Crumley E, Wiebe N. Tissue adhesives for traumatic lacerations in children and adults. *Cochrane Database Syst Rev.* 2002;(3):CD003326. Review.
8. Haghpanah S, Vafafar A, Golzadeh MH, Ardeshiri R, Karimi M. Use of Glubran 2 and Glubran tissue skin adhesive in patients with hereditary bleeding disorders undergoing circumcision and dental extraction. *Ann Hematol.* 2011 Apr;90(4):463-8
9. Amiel GE, Sukhotnik I, Kawar B, Siplovich L. Use of N-butyl-2-cyanoacrylate in elective surgical incisions--longterm outcomes. *J Am Coll Surg.* 1999 Jul;189(1):21-5.
10. Topical Skin Adhesive: Strong, Safe, Secure Wound Closure. Ed. Pfiedler Enterprises 2101 S. Blackhawk St, Suite 220, Aurora, CO 80014 www.pfiedlerenterprises.com
11. Hosseini SM, Rasekhi AR, Zarenezhad M, Hedjazi A. Cyanoacrylate glue dressing for hypospadias surgery. *N Am J Med Sci.* 2012 Jul;4(7):320-2.
12. Howell JM, Bresnahan KA, Stair TO, Dhindsa HS, Edwards BA. Comparison of effects of suture and cyanoacrylate tissue adhesive on bacterial counts in contaminated lacerations. *Antimicrob Agents Chemother.* 1995 Feb;39(2):559-60
13. Amin M, Glynn F, Timon C. Randomized trial of tissue adhesive vs staples in thyroidectomy integrating patient satisfaction and Manchester score. *Otolaryngol Head Neck Surg.* 2009 May;140(5):703-8.
14. Data on GEM S.r.l. file.
15. Singer AJ, Quinn JV, Hollander JE. The cyanoacrylate topical skin adhesives. *Am J Emerg Med.* 2008 May;26(4):490-6.
16. Osmond MH, Klassen TP, Quinn JV. Economic comparison of a tissue adhesive and suturing in the repair of pediatric facial lacerations. *J Pediatr.* 1995
17. S Masi, L Mugnaini, A Pazzaglia, C Saffirio, A Salvatore, I Sforzi, F Silenzi and S Gianassi. Synthetic Two-Components Skin Glue (Glubran Tiss) Effectiveness and Suitability in Pediatric Emergency Department. *Pediatric Emergency Care and Medicine: Open Access.* 2017 - Vol. 2 No. 1:1.
18. Resistenza alla trazione misurata al nodo. The United States pharmacopeia. Suture. USP 41: United States Pharmacopeial Convention; 2017. p. 3902.

A range of solutions in continuous evolution.



➤ G-NBOC2-25



➤ G-NBOC2-35



➤ G-NBOC2



10 sterile single doses in aluminium blister packs - 2 year validity

**SOLUTION
COMES FROM
EVOLUTION.**



Via dei Campi 2 - PO Box 427 - 55049 Viareggio (LU) Italy
Tel. +39 0584 389784/391388 - Fax +39 0584 397904
www.gemitaly.it - info@gemitaly.it